Table of Contents

DRF™ (Dynamic Radio Frequency) Technology – How does it work?	4
Intended Use	5
Who Should NOT Use the Device?	6
CONTRAINDICATIONS	6
WARNINGS	7
CAUTIONS	8
Possible Adverse Effects	9
Device Description	10
User Instructions	12
Preliminary Test	12
Getting Started	13
Recommended Treatment	14
Cleaning & Storage	17
Troubleshooting	17
Customer Service	18
Device Specification Parameters	19
Labelling	20

Introducing DRF™ (Dynamic Radio Frequency) Technology – How does it work?

This device brings a new approach to non-invasive skin rejuvenation featuring Sensica proprietary Dynamic RF™ technology, combining clinically proven RF (Radio Frequency), commonly used by doctors, with advanced facial massage.

Radio Frequency (RF) is, a form of electromagnetic energy, intended to reduce wrinkles, tighten skin and renew skin elasticity. It is suitable for all skin tones, easy to use, and provides noticeable, long-lasting results.

RF technology has been used in doctors' clinics and spa centers for several years and was found to be highly effective. The device has two electrodes that deliver the RF energy into the skin to contract collagen fibres and stimulate the natural production of collagen, which decreases as we age. If used correctly and according to the

Recommended Treatment Protocol, the device will reduce the appearance of wrinkles as well as tighten skin and renew elasticity. Results may be visible after just one treatment.

Please make sure to carefully read this user guide before you start treating with the device.

The user guide can also be found at www.sensica.com

Intended Use

The device is intended:

- For women and men who want to: reduce facial wrinkles, tighten their facial skin, and renew skin elasticity.
- For use on: face, cheeks, "smile" folds, "crow's feet", forehead, chin, jawline, allowed areas of the neck, décolletage, and hands.
- For use on fine to moderate depth wrinkles.
- To be used with the Base Gel. This gel is especially designed to work with Sensica RF devices and ensures proper delivery of RF energy. This will provide the best possible results by improving contact between the device and the skin to improve skin conductivity.

Note: It is required to apply the Base Gel in order to correctly use the device. If the Base Gel is not applied, it will compromise your safety and treatment efficacy.

Who Should NOT Use the Device?

Carefully read the information below before you begin your first treatment. If you still have any doubts or are unsure that the device is safe for your use, please consult with your physician.

CONTRAINDICATIONS

Do NOT use the device if you:

- Are under the age of 18. The use of Sensilift on or by children has not been established and could result in serious injury.
- Have an active implant such as a pacemaker or insulin pumps or if you have any severe known heart disorder such as arrhythmias.
- Are pregnant or breastfeeding.
- Suffer from severe or uncontrolled endocrine or metabolic conditions (such as diabetes and uncontrolled thyroid disease).
- Have a damaged immune system or autoimmune disorder such as AIDS, or Lupus, or if you use drugs that influence the immune system.
- Suffer from blood clotting disorders or excessive bleeding, bruising or damaged wound healing or if you are routinely taking blood thinning medications.
- Suffer from a history of deep vein thrombosis (formation of blood clots in the veins).
- Have a history of collagen disorders or suffer from keloid scar formation (growth of extra scar tissue where the skin has healed after an injury).
- Have an extreme general illness, such as hepatitis, liver disease or epilepsy.
- Suffer from sensitive skin such as in cases of: frequent redness, inflammations, breakouts, infection, acne or extra dry skin.
- Have diminished or exaggerated perception of temperature changes.
- Have been taking Accutane in the last 6 months.
- Have a heat induced disease or an active skin disease such as: Herpes, Eczema, rash etc. in the treatment area.
- Have skin cancer, and do not use on areas at risk of malignancy or if you have received radiation therapy or chemotherapy within the past 3 months.

WARNINGS

Do **NOT** use the device:

- Over any permanent implants in the designated treatment area including, but not limited to: silicon or injected chemical materials, metal plates or other metal materials.
- If you had plastic surgery, such as face lifting or eyelid surgery, or if you had aesthetic injections of fillers such as: Hyaluronic Acid, Collagen, fat or Botulinum Toxin, in the treatment area within the past 3 months.
- If you have had a professional chemical peel within the past 3 months with Alpha Hydroxyl Acids (AHA), such as glycolic acid, or Beta Hydroxyl Acid (BHA), or if you had a professional facial abrasion such as: dermabrasion, facial resurfacing, deep chemical peeling (phenol, Trichloroacetic acid-TCA), and/or ablative laser treatments within the past 3 months.
- Over nostrils, lips of the mouth, ears, breast or genital area, or areas not recommended in this user manual.
- Over the thyroid since it may cause adverse events.
Over the eyelids, eyes or sides of the eyes (next to your nose), since it may cause eye damage.
- Over any metal items (such as earrings, jewelry or piercing) or on tattoos or permanent make up in the designated treatment areas.
- Over sunburns, dark brown or black spots such as: freckles, birth marks, around the nipples, moles or blisters in the treatment area.
- If you are allergic to any components of the Base Gel.
- If you have a sensitivity to chrome or nickel as these components are found in the two electrodes of the device.

Additionally, this device is not intended for use:

- By persons with compromised physical, sensory or mental capabilities.
- Without removing excessive facial hair from treatment areas, since it may result in burns in these areas.

Note: If you have undergone any other aesthetic procedures in the designated treatment area or have any uncertainty regarding potential side effects, consult your physician before using the device.

CAUTIONS

• KEEP OUT OF REACH OF CHILDREN.

Never use flammable liquids such as alcohol or acetone, including items like perfumes or sanitizers, on the skin one hour before using the device as it may cause burning of the skin.

- Do not open or repair the device as it may expose you to dangerous electrical components.
- Do not place, store or use the device near: bathtubs, showers, basins or other vessels containing water as this may cause electrocution.
- Do not leave the device unattended while plugged into the electricity as it may cause electrocution.
- Use the device only according to the Instructions

- Do not use the device if the device is not working properly or any part of the device is missing, damaged or loose. Stop use if you see or smell smoke or if the device becomes damp or wet.
- Do not use the device without applying the Base Gel before use.
- Do not use the Base Gel after its expiration date.
- Before using the device, make sure it is clean, dry and intact. Cleaning instructions are detailed in the cleaning section.

Note: this device is designed to be used by a single user.

Possible Adverse Effects


When using the device according to the Instructions for Use, adverse effects are rare. However, every aesthetic procedure, including the usage of home-use devices, may cause some degree of adverse effects (as listed in the table below).

In case adverse effects appear, stop using the device and consult with your physician if necessary.

Mild adverse effects should be spontaneously relieved within a few hours.

Severity of Adverse Effect	Frequency of Adverse Effect	Adverse Effects
Mild	Common	Local pain Skin irritation Skin redness Skin sensitivity Allergic contact dermatitis or sensitivity to the Base Gel
Moderate	Rare	Mild skin burns Pigment changes Excessive redness and swelling Urticaria, purpura or ecchymosis Hematoma
Severe	Extremely Rare	Severe burns Infection Bruising Blistering

Device Description

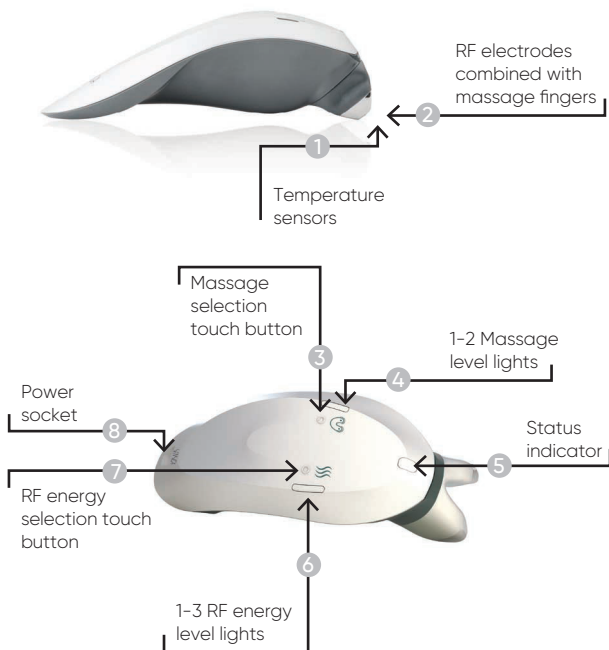User interface description

1. ON/OFF and energy selection button – A gentle touch on the button (after the device is connected to the power supply) will turn the device ON at energy level 1, shown by 1 illuminated light in the energy level display. Once it reaches the highest level, another gentle touch will turn the device OFF. The device will switch to standby mode when left unattended for 1 minute.

2. Status indicators

Status Indicator	Blinking yellow light	The device is on Stand-by mode . In case of no light, make sure the power cord is plugged-in properly.
	Steady green light	The device is in Treatment mode . RF is activated, skin contact is detected and the treatment is undergoing.
	Blinking alternating green and yellow lights	The device is in Error mode . the device detected a malfunction and prevents the activation of the RF and massage functions.
RF energy level selection touch button 	1-3 energy levels lights	Green light on, indication of RF energy level. Touch the RF energy level selection button repeatedly until you reach the desired level: 1, 2, 3 or off.
Massage level selection touch button 	1-2 massage level lights	Green light on, indication of massager function level. Touch the massage level selection button repeatedly until you reach the desired level: 1, 2, or off.

User instructions

The device is intended for use on: cheeks, nasolabial "smile" folds, "crow's feet" (next to your eyes), forehead, chin, jawline, allowed areas of the neck, décolletage and hands.

What to Expect

During your facial wrinkle reduction and skin tightening treatment it is normal to experience:

- Sensations of warmth during the treatment session – created by the RF energy during treatment.
- Some mild red or pink color on the skin – during and just after each treatment.

Preliminary Test

The preliminary test is intended to determine your skin's reaction to the treatment and help you find the best settings for each treatment area. Perform a preliminary test on each area you intend to treat 48 hours prior to the full treatment.

- Prepare your skin as described below in "**Step 1: Prepare your skin**". Instructions for operating the device can be found below in "**Step 2: Device operation**".
- Before using the device, apply an even layer of the Base Gel on the treatment area.
- Start your test with energy level 1 and massage level 1, and move the device in small circular motions.
- Increase the levels according to your comfort. Stop the test if you do not feel comfortable or if you see adverse effects.
- Perform the preliminary test for a total time of about 3 minutes.
- Wait 48 hours and examine the tested area.

- If the skin in the tested area appears normal and no adverse effects are present, you can proceed with a full treatment at the highest comfortable energy level that you tested.

Note: The RF energy level indicates the rate at which the temperature rises when in heating mode. A high RF energy level will enable faster and more even heating. You can achieve maximum results with all RF energy levels but for lower levels it will take more time to reach the effective temperature. Therefore, you may select any RF energy level that you find comfortable.

Getting Started

Step 1: Prepare your skin

1. Remove jewelery from the treatment area, as well as from your hands, prior to use.
2. Prior to use, make sure your skin is dry, clean and free of topical products such as: moisturizers, liquids or lint.
3. Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin one hour before using

the device as it may cause burning of the skin.

4. The device is suitable for men. However male users should make sure their skin is clean-shaven before treatment.

Step 2: Device Operation

Note: A preliminary test should be conducted separately for each area you intend to treat before proceeding with a full treatment in that area. See **"Preliminary Test"** section above.

1. Make sure the device is clean, dry and intact prior to use.
2. Apply an even layer of "Base Gel", included in your kit, on the area of your skin intended for treatment.
3. Set the Radio Frequency (RF) energy level and the massage speed to the level you selected in your preliminary test treatment. Gently touch the RF energy level selection button  repeatedly until you reach the desired level: 1, 2, or 3. Gently touch the massage level button  repeatedly until you reach the desired level: 1 or 2.

4. The activation of either RF or massage functions can be done independently or together. Both functions will be activated at level 1. Select different levels based on the preliminary test.
5. Sensilift is equipped with 2 safety sensors which will keep the temperature in the effective and safe range during the whole treatment.
6. Place the treatment fingers on the intended area, ensuring full contact with your skin during the entire treatment. Sensilift will indicate the proper skin contact by the green status light. If Sensilift is not properly placed, the RF energy will not be delivered.
7. If you do not feel sensations of warmth after approximately 10 seconds, try moving Sensilift on a smaller treatment area, and/or increasing the RF energy.
8. Move the device on the treated area in slow circular or back and forth motion while constantly and firmly contacting the skin.
9. See treatment areas, treatment time and treatment schedule in the "Recommended Treatment Protocol" chapter.
10. After completing the treatment, wipe off the remaining "Gel Base" from the treated area.
11. Cleaning instructions for before and after are depicted in the Cleaning chapter.

Recommended Treatment

For long-lasting results, use the device once a week on each desired area for a

- period of 8 weeks.

Response to the treatment differs between individuals and depends on

- environmental and biological factors, heredity and aging.

To sustain optimal long-term results, maintenance treatments may be

- performed as needed every 4-8 weeks.

Treatment areas

Do not treat full areas with energy levels you have not tested before

Forehead (treatment area 1)

- Place the treatment fingers on the area as indicated in the picture.
- Keep moving the device up and down across the forehead area.
- Recommended treatment time for this area is 10 minutes.

Between the eyes – "Eleven" (treatment area 2)

- Place the treatment fingers on the area as indicated in the picture.
- Keep moving the device in small circular movements or up and down
- Recommended treatment time for this area is 10 minutes.

"Crow's feet" (around the eyes) (treatment area 3)

- Place the treatment fingers on the area as indicated in the picture.
- Keep moving the device in continuous, slow 'side to side' circular motions.
- Avoid contact with the eyes and eye cavities, as treatment in these areas may cause adverse effects.
- Recommended treatment time for this area is 5 minutes (for each side).

Cheek and nasolabial "smile" folds (treatment area 4)

- Place the treatment fingers on the area as indicated in the picture.
- Move the device in slow circular motions across the cheek area up and down over the "smile folds".
- Avoid areas such as inside the nostrils and oral cavity, as treatment in these areas may cause adverse effects.
- Recommended treatment time for this area is 5 minutes (for each side).

Chin treatment (treatment area 5)

- Place the treatment fingers on the area as indicated in the picture.
- Keep moving the device from one side of the chin to the other.
- Recommended treatment time for this area is 5 minutes.

Neck treatment (treatment area 7)

- Place the treatment fingers on the area as indicated in the picture.
- Keep moving the device up and down across the neck area.
- The recommended treatment time for this area is 5 minutes (for each side).

Caution: Avoid treatment in the lower third of the neck since this area is more sensitive, avoid from treating the thyroid area (marked .

Jawline treatment (treatment area 6)

- Place the treatment fingers on the area as indicated in the picture.
- Keep moving the device up and down along the jaw.
- Recommended treatment time for this area is 5 minutes. (for each side).

Décolletage and hands treatment

- Place the treatment fingers on the desired treatment area.
- Keep moving the device in slow circular motions across the décolletage or hands area.
- Since décolletage and hands are large areas, for best results it is recommended to divide those areas to sub-areas of approximately 4cm².
- Recommended treatment time for each sub-area is 5 minutes.

Post facial treatment care

- We recommend you also use Sensica's anti-aging cream to increase treatment efficacy and ensure long-lasting results (available for purchase at www.sensica.com).
- Apply a high-factor sunscreen (at least 30 SPF) and protect the treated area from sunlight for two weeks following treatment.
- Following the treatment, avoid any usage of hot water on your face and discontinue usage of abrasive or harsh products for two days.
- Makeup may be applied immediately after treatment.

Cleaning & Storage

- Make sure to turn off and unplug the device before cleaning.
- It is recommended to clean the device, especially the treatment fingers, temperature sensors, and silicone cover, of residual from the Base Gel with a soft damp cloth or paper towel.

Do not rinse the device under running water, as it may cause device malfunction.

- Keep the device dry. Never place, store or use the device near bathtubs, showers, basins or other vessels containing water as this may cause severe electrocution.

Troubleshooting

If you experience any problems with your Sensilift, please refer to the information below to try to resolve it:

- If the device does not turn on:
 - a. Make sure that the power cord is plugged properly into the electrical outlet and into the device power socket.
 - b. Make sure the power outlet has power.
- If the device does not deliver RF energy, meaning that you cannot feel the warm sensation, be sure to:

Verify there is full contact with the skin and that the Base Gel was applied properly on the skin. The RF

signal will be activated only when the treatment fingers are in full contact with the skin. See indicator states in User Interface description chapter.

- In case of dirt or debris stuck on the treatment fingers, you should clean them in accordance to the Cleaning chapter. If you can't get the area clean, contact our customer service.
- In case of a general error, reset the device by unplugging the power cord from the electrical outlet on the wall, wait 10 seconds and then reconnect it. If you still have problems with your device, please contact our customer service.

Customer Service

If your device has broken/stopped working or if you need any further assistance/information, please contact us. You may also contact us to order additional Base Gel or you can purchase it online at sensica.com.

 **+44-800-923-4830**

 **supportUK@sensica.com**

Device specification parameters

The RF handheld device is powered by bipolar RF at 10 ± 0.5 [Watt]

Frequency: 1 ± 0.05 [MHz]

Selectable output power level:
3 levels

Mechanical massage:
implemented by moving RF electrodes. Speed: 50-100 RPM

Temperature sensors: RF heating is automatically deactivated when skin temperature reaches 40°C (104°F). The heating is reactivated when temperature decreases.

Treatment area: ~ 2.2 [cm^2] / 0.34 [in^2]

External power supply: Input 100-240 [VAC], 50-60 [Hz], Output 9[VDC], 1.5[A]

Device Size and Weight:
140 [g], 173[mm] X 56.8[mm] X 53.1[mm] LxHxW
4.93 [Oz], 6.8[in] X 2.23[in] X 2.09[in]

Storage Instructions

Storage and Transport Conditions

Temperature	-25°C to 70°C (-13°F to 158°F)
Humidity	10% to 90% RH at $55^{\circ}\text{C}/131^{\circ}\text{F}$
Pressure	Non-condensing

Operational Conditions

Temperature	5°C to 40°C (41°F to 104°F)
Humidity	20% to 80% RH
Max Operating Altitude	3,000 m.

Labeling

Model: STG100EU

Power: 9 [V], 1.5 [A], 100–240 [VAC], 50–60 [Hz]

Triple insulation

Consult Instructions for use

Waste Electronic and Electrical Equipment

Please read all warnings and safety precautions prior to use

Do not rinse the device under running water, as it may cause device malfunction. Keep dry

Use the device only with original power supply. For your safety, never use DC adapters which are not approved by the manufacturer.

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