Table of Contents

RPL™ (Reactive Pulsed Light) Technology-How does it work?	2
Intended Use	3
Who should NOT use the device?	4
CONTRAINDICATIONS	4
WARNINGS	5
CAUTIONS	6
Possible Adverse Effects	7
Device Description	8
User interface description	9
User Instructions	12
Preliminary Test	12
Getting Started	13
Step 1: Prepare your skin	13
Step 2: Device operation	13
Cleaning & Storage	19
Troubleshooting	19
Customer Service	20
Device specification parameters	21
Labeling	22

RPL™ (Reactive Pulsed Light) Technology - How does it work?

RPL™ (Reactive Pulsed Light) is an advanced technology which uses clinically proven IPL (Intense Pulsed Light) technology in conjunction with a skin tone and skin contact sensor. The light emitted from the device penetrates the skin and is then absorbed by the melanin pigment in the hair shaft. The heat produced by this light goes into the follicle and generates a local thermal effect, also known as selective photothermolysis. The thermal effect disables the hair follicle which causes a reduction in hair growth while leaving the surrounding skin unaffected and cool.


Long-term hair reduction- Hair growth cycle

Two factors play an important role in optimal long-term hair reduction:

- **Hair growth cycle.** Hair growth involves three phases: Anagen, Catagen, and Telogen. Anagen is the

growth phase, while Catagen and Telogen are the resting phases. Permanent hair reduction can only occur in the Anagen phase, the phase of growth. Since hair growth is not synchronized, it is recommended to complete the full treatment protocol to ensure every hair will be exposed to treatment during its Anagen phase.

- **Skin tones (according to the Fitzpatrick scale).** Darker skin contains more melanin, which competes with the targeted hair for energy absorption. When more energy is absorbed into the skin, there is a greater risk of adverse effects.

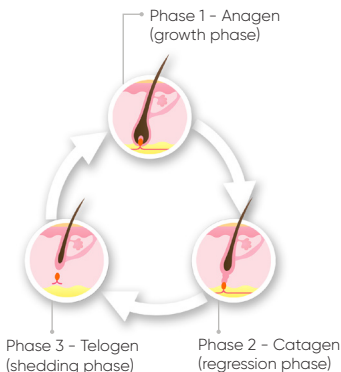Please make sure to carefully read this user guide before you start using the device. The user guide can also be found at www.sensica.com

Intended Use

This is a home-use device intended for the removal of unwanted hair on the entire body (e.g., legs, arms, underarms, bikini line, back, shoulders, chest and stomach) and on permitted facial areas (i.e., above upper lip, cheek, jaw, chin and neck).

It is considered a medical device without medical use and does not have clinical benefit.

Note: This device is especially designed for your own personal treatment and is intended for a single user and home use only.

The device is intended for use on skin tone types I to V according to the Fitzpatrick skin chart below (coloured chart available on the cover). Do not use on dark skin such as type VI. The use

of IPL devices on naturally dark skin is forbidden.

Skin Tones

Hair Colour

This device works best on dark hair or hair that contains more melanin. Brown and light brown hair will respond, but typically will require more treatment sessions. Red hair may show some response while white, grey or blond hair will show a limited response. Nevertheless, some users have noted hair reduction after multiple treatment sessions.

Note: Hair removal by lasers or Intense Pulsed Light can cause increased hair growth in some individuals. Based upon current data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage, treated on the face and neck.

Who should NOT use the device?

Carefully read the information below before you begin a treatment with the device. If you still have any doubts or are unsure that the device is safe for your use, please consult with your physician.

CONTRAINDICATIONS

Do NOT use the device if you:

- Have an active implant such as a pacemaker, insulin pump or similar devices, since it may cause interference.
- Are pregnant, have a chance of being pregnant or are nursing since the device was not tested on pregnant or nursing women and may stimulate pigmentation.
- Have skin cancer, areas at risk of malignancy, or if you have received chemotherapy/radiation therapy within the past 3 months since the safety of such usage has not been established and the immune system can be compromised, which can lead to serious injury.
- Suffer from epilepsy, as using the device could trigger a seizure.
- Have damaged or broken skin, or a history of Herpes or Psoriasis outbreak in the treatment area since such usage can result in serious side effects.
- Suffer from keloid scar formation (growth of extra scar tissue where the skin has healed after an injury) in the treatment area since such usage has not been established and may result in serious skin complications.
- Have been taking Tetracycline antibiotics or medication with the following active ingredients within the last two weeks: Retinoid acid, Retina-A, Vitamin A (such as Accutane), as they are considered to be photo-sensitizing drugs and may result in skin injury.
- Are photosensitive or have other diseases related to light sensitivity, as device usage under these conditions may lead to serious injury.

- Suffer from severe or decompensated metabolic conditions, such as diabetes, since it may cause a greater risk of delayed healing and infection in case of an adverse effect.
- Have had a skin peel or other skin resurfacing procedures within the last 8 weeks in the treatment area as, under these circumstances, device usage may result in stimulation of pigmentation.

Note: If you have undergone any other aesthetic procedures or have any uncertainty regarding potential side effects, consult with your physician before using the device.

WARNINGS

Do NOT use the device:

- On naturally dark skin (coloured chart available on the cover).
- If you have dark brown or black spots such as: around the nipples, freckles, birth marks, moles, or blisters on the area to be treated.
- On mucous membranes such as: nostrils, lips of the mouth, eyelids, ears, genital area and the anus, as it may lead to serious injury.
- Over piercings or any metal items such as earrings or jewelry on the area to be treated, as it may cause burns. Over areas with tattoos or permanent make-up, cover tattoos or permanent make-up with a white plaster before treatment. Do not use the device directly over the tattoo or permanent make-up, even when covered.
- **On tanned skin or after sun exposure as this can cause serious burns or skin injury. Avoid more than 15 minutes of sun exposure for 4 weeks before your treatment and for at least 2 weeks post treatment. If during this period your treated skin area is exposed to sunlight, make sure to apply SPF 30 or higher sunscreens.**

Additionally:

- Never look directly at the light coming from the applicator and lamp, as it may cause eye injury.
- Never use flammable liquids such as alcohol (including

- perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin one hour before using the device as it may cause burning of the skin.
- Do not treat the same area of skin more than once per treatment, as it may cause adverse effects.
 - Do not leave the device unattended while plugged into the electricity, as it is a risk of electrocution.
 - Do not use the device if any part of the device or its accessories are: not working properly, appear damaged, broken, loose, missing or wet (take extra care since UV and IR are emitted from this device), or if you see or smell smoke.

CAUTIONS

- KEEP OUT OF REACH OF CHILDREN.
- This device is not intended for use by individuals with reduced physical, sensory or mental capabilities.
- Do not place, store or use this device near bathtubs, showers, basins or other vessels containing water, as it may lead to severe electrocution.
- Do not open or repair the device, as it may expose you to dangerous electrical components which may cause serious bodily damage.
- Before using this device, make sure it is clean, dry and intact. Cleaning instructions are detailed under the **Cleaning & Storage** section.

Possible Adverse Effects

When using this device according to the Instructions for Use, adverse effects are rare. However, every aesthetic procedure may cause some degree of adverse effects (as listed in the table below). In case adverse effects appear, stop using the device, and consult with your physician if necessary.


Mild adverse effects should be spontaneously relieved within a few hours.

Severity of Adverse Effect	Frequency of Adverse Effect	Adverse Effects
Mild	Common	Skin Irritation Skin Redness Skin Sensitivity
Moderate	Rare	Skin Burns Scarring Pigment Changes Excessive Redness and Swelling
Severe	Extremely Rare	Infection Bruising

Any serious incident that has occurred in relation to the device should be reported to us (via customer service) and to the competent authority of the Member State in which you are a resident.

Device Description

1. Flashing light (internal lamp)
2. Treatment window
3. Trigger button
4. Energy level display
5. Battery status indicator
6. Energy level selection buttons
7. Energy recommendation button
8. ON/OFF button
9. Status Indicator
10. Skin tone and skin contact sensor
11. Charging socket
12. Large treatment window (0.7 sq. in/ 4.5 cm²)
13. Precise treatment window (0.3 sq. in/ 2 cm²)
14. Power adapter

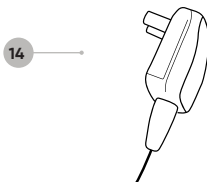User interface description

1. **ON/OFF button** – Pressing the button will turn the device ON at energy level 1. The device will switch to standby mode when left unattended for 5 minutes or by pressing the ON/OFF button for 3 consecutive seconds.
2. **Energy recommendation button** – A press on the **Q** button will allow you to detect the recommended energy level for your treatment, according to your skin tone in the tested area. When pressing the skin tone detection button, it will light up green and all energy level display indicators will flicker. Hold the device on the desired treatment area and press the trigger button to measure your skin tone (no flash will be released). The recommended energy level according to your skin tone at the treatment area will be automatically set, as indicated by the energy level display. Failure to determine suitable energy level due to lack of contact or not allowed skin tone

(too dark), is indicated by a momentary red light on the status indicator. Note that the recommended energy level is set only through this procedure and it does not occur automatically. To test the skin tone in each treatment area, place the device on that area and press the energy recommendation button. Performing a preliminary test in the desired treatment area is still needed and selection of energy level should be according to your own comfort, as determined by the preliminary test.

3. **Energy level selection buttons** – The + (increase energy level) and – (decrease energy level) buttons, allow you to choose between 6 different levels. Select the energy level according to what you found most comfortable in the preliminary test. You can use the energy recommendation button to assist in choosing the optimal energy level for you for each desired treatment area.

- Each press of the (+) or (-) buttons will increase or decrease the energy by one level. When you reach the maximum energy level, pressing the (+) button again won't change the energy level any further. Similarly, when you reach the minimal energy level, pressing the (-) button won't decrease the energy level further.
- The 1-6 energy level display indicators will light up green according to the selected energy level (e.g. at energy level 1, one green light will be lit, at energy level 2, two green lights will be lit).

4. Battery status indicator –

The indicator, marked as , shows the charge status of the battery:

- Red – Battery is being charged. Note: If in addition to this red light, energy level "1" is green, the device is turned on.
- Green – Battery is fully charged, and the device is connected to the external power supply.
- Slow blinking red light – Low battery level (about 30% of full charge).
- Fast blinking red light – Very low battery level (about 10% of full charge).

Note: The green or blinking red battery status indicators will appear only while the device is turned on.

5. Treatment windows – For the most comfortable treatment, select the appropriate treatment window for the body area you want to treat and simply snap on the treatment window.

- **Large treatment window (the optical filter is fixed on the device)** – Use the large treatment window for fast and easy hair removal on large body areas, such as legs, back and arms.
- **Precise treatment window (additional snap-on window)** – Snap on the precise treatment window for treatment of small and delicate areas, such as the face (area below the cheek line) and bikini line.

6. Status indicator

Status Indicator	Status	Description
No light	OFF	Device is inactive
Steady green light	ON	Flash is ready to be delivered
Momentary orange light	Flash blocked	The device is not in full contact with the skin – no flash will be delivered
Momentary red light	Flash blocked	A non-allowed skin tone has been detected – no flash will be delivered (coloured chart available in section 02 on the cover)
Blinking orange light	Energy recommendation	Device ready for skin tone measurement – no flash will be delivered during this time
Steady red light	Device blocked	The device is blocked due to either an internal malfunction or improper attachment of treatment window to the device

User Instructions

The device is suitable for use on: legs, arms, underarms bikini line, back, shoulders, chest, stomach, and permitted facial areas (above upper lip, cheek, jaw, chin and neck).

What to expect

During your treatment, it is normal to experience:

- A flash of light with each press on the trigger button– the emitted light will not harm your eyes.
- A pop noise with each flash of light.
- A mild, stinging sensation.

Preliminary Test

The preliminary test is intended to determine your skin's reaction to treatment with the device and help you set the optimal energy level for each treatment area.

Perform a preliminary test on a hidden patch of skin for the area you intend to treat 48 hours prior to the first full treatment.

Due to varying levels of sensitivity, body areas may require different energy levels; this will be determined by separate preliminary tests.

Note: Do not treat full areas with energy levels that you have not tested before!

Do not treat the tested area for at least 48 hours after the preliminary test has been carried out.

1. Prepare your skin as described below in "**Step 1: Prepare your skin**". Instructions for operating the device can be found below in "**Step 2: Device operation**".
2. Select the lowest energy level and deliver a flash.
3. You can increase the energy level as long as it feels comfortable. You may continue to do this up to the maximum energy level unless you begin to experience discomfort.
4. If you do not feel comfortable or see adverse effects such as severe redness or blistering stop the test immediately.

5. Wait 48 hours and examine the tested area:
 - If the skin in the test area seems normal (there is no reaction or just a mild redness), you can proceed with a full treatment on the highest comfortable energy level that you tested.
3. Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone on the skin for one hour before using the device, as it may cause burning of the skin.

Note: You may use the Energy Recommendation button to automatically detect the recommended energy level for each treatment area. Please note this is a recommendation based on your skin tone and you are recommended to test that this energy is comfortable for you.

Getting Started

Step 1: Prepare your skin

1. Remove all jewelry.
2. Make sure the hair in the treatment area is shaved or trimmed (between treatments you can shave, but don't wax or pluck hair), and that the skin is free of any residue from: powders, creams, makeup, antiperspirants or deodorants.

Important: Do not use any cream or gel while using the device since it may prevent normal operation.

Remember: The device is equipped with a built-in skin-tone sensor, allowing the device to work only on suitable skin tones. This unique safety feature will prevent a flash from occurring on specific areas where your skin is too dark or too tanned, even if your general skin type is lighter.

Step 2: Device operation

Note: A preliminary test should be conducted before proceeding with a full treatment. See "Preliminary Test" section above. Make sure the device is clean, dry and intact. You can use the device while it is being charged. However, this may limit the number of pulses since battery charging is slower than battery consumption during use.

1. Charge your device as necessary:
 - Before your first use and when the battery is empty, fully charge the battery (full charging of the battery takes up to 2.5 hours).
 - Whenever battery level is low, as indicated by the blinking red light of the battery status indicator.
 - Fully charge the appliance at least once every 6 months, even if you have not used it in that time.
 - To charge your device, insert the plug of the power cord into the electrical outlet and the other end into the device. When the battery is fully charged, the charging light lights up green continuously.
 - After charging, remove the power adapter from the wall socket and disconnect it from your device.

Notes about charging:

- Some warming of the device and power adapter during charging is normal.
 - Charging should be done in a room temperature or cool environment, never above 40°C. The device is equipped with battery-overheat protection and will not charge if the room temperature exceeds 40°C.
 - Never cover the appliance and the adapter during charging.
2. Turn the device on. When the device is turned on it will be automatically set to the lowest energy level and the energy level display will show one green light. You may use the energy recommendation button at the start of each treatment or manually change energy levels by pressing the energy level selection buttons (increase energy level button (+) or decrease energy level button (-)). You may decrease or increase the energy level according to your comfort level, as determined by the preliminary test.

3. Select the desired treatment window and simply snap it onto the device by aligning the plastic pins of the selected window to the device sockets and pressing into place. To remove the treatment window, grab both sides on the elevated dots and firmly pull until it is removed.
4. Place the treatment window on the desired treatment area, ensuring full contact with the skin. As a safety precaution the device is equipped with a safety sensor that will not trigger a flash if:
 - The device is not properly applied to the skin. In this case, the status indicator will turn "momentarily" orange.
 - Your skin tone is too dark/tanned. In this case, the status indicator will turn "momentarily" red.
 - The treatment window is not properly attached. In this case, the status indicator will turn "steady" red.
5. Apply the treatment window firmly against your skin and press the trigger button to release a flash. Following each flash, move the device to the next spot. After each flash, the device will take 1-3 seconds to be ready for the next one.
6. You may choose from 2 operational modes, the 'Pulse mode' or the 'Glide mode':
 - **Pulse mode:** Press the trigger button once to release a single flash. Continue to the next spot and press the trigger button again to release another flash. Repeat this until treatment is complete.
 - **Glide mode:** Begin treatment in Pulse mode. After three consecutive flashes without releasing the trigger button, the device will automatically enter the Glide mode. This mode will release flashes continuously, with a 1-3 second pause between each one, so there is no need to hold the trigger button. The Glide mode will automatically

deactivate when the device is not in proper contact with the skin or if the skin tone is too dark.

7. After the treatment, turn off the device by pressing the "ON/OFF" button; the status indicator will turn orange while the device shuts down.
8. Perform treatment sessions according to the **"Recommended Treatment"** section.

Important note: Try not to overlap pulses – repeated pulses on the same spot may cause burns.

Recommended Treatment:

Body hair:

Hair removal treatments **1–4** will be two weeks apart.

Hair removal treatments **5–8** will be four weeks apart.

Facial hair:

Hair removal treatments **1–5** will be two weeks apart.

Hair removal treatments **6–8** will be four weeks apart. To maintain the results following completion of the treatment

protocol, we recommend performing a treatment once a month, or as required.

What to expect:

Different individuals and different treatment areas may not react equally to RPL technology.

You should begin to notice a reduction in hair growth by the 6th treatment. Once the treatment protocol of 12 treatments is complete, you should see further hair growth reduction. Any hairs that remain should be finer and lighter in colour. In order to maintain results, we recommend using the device once a month per treatment area.

Important: Facial treatments should be performed at the medium or low energy level. Don't forget to select the proper treatment window (either the built-in large treatment window or the additional precise treatment window), according to the size of the treatment area and your comfort.

Notes:

- It is recommended to use the highest energy setting that you found comfortable during the preliminary test.
- Lighter hair colors may require more treatments in order to see results.
- The hair will continue to grow for up to two weeks after a hair reduction session. However, as you progress through the full treatment protocol, hair will naturally and gradually fall out. The time until initial results can be noticed varies between users and depends on factors such as: age, hair color and skin tone.
- The number of sessions required for achieving long term results differs for every treatment area and between different users.
- This device is suitable for men. However, male body hair, especially on the chest, will require more treatment sessions.

- Maintenance will be needed from time to time if growth is still visible.

Facial treatment

The facial skin is a sensitive area that requires special care and caution. Make sure you follow all instructions and safety precautions carefully in order to avoid adverse effects. It is recommended you first do a body treatment and get to know the device before performing facial treatments.

This device can be used to treat all facial areas below the cheek line, such as: above upper lip, cheek, jaw, chin and neck.


Never use the device around or near your eyes.

Warnings!

- When using this device to treat facial areas, do not use energy levels exceeding the medium level, since it may cause severe skin damage.
- Do not use the device on the face above the cheek line or around the eyes, eyebrows and eyelashes as this could cause serious eye damage.
- Avoid the nostrils and lip, as these areas are more sensitive.
- The skin on the lower third of the neck is more sensitive. Take extra care to avoid overlapping treatment areas.

Note: Facial hair may be hormonal and as such, it can take more sessions to achieve a permanent hair reduction.

For male users: Take care while treating the beard area; the density of hair means a large amount of melanin which absorbs the heat during treatment.

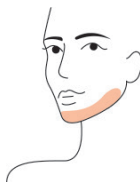Post-treatment care

- After treatment, you may use a calming lotion.
- Avoid using peeling and bleaching creams or similar products for at least 24 hours, as this irritates the treated skin.
- Avoid direct sun exposure for at least 2 weeks post treatment. If exposing recently treated skin to sunlight during this time period, apply sunscreen of SPF 30 or higher.

Cleaning & Storage

1. Make sure to turn off the device before cleaning and storage.
2. It is recommended that the treatment window area be cleaned after every session by using a dry and clean cloth.
3. Fully charge the device every 6 months, even if you have not used it in that time.
4. If there is dirt or debris on the treatment window,

clean the treatment window. If you can't get the device clean, contact our Customer Service.

Troubleshooting

If the device does not turn on:

- Make sure that the device is charged (as indicated by the green light at the battery status indicator).

If the device does not trigger a pulse:

- Verify there is full contact between the device and your skin. The trigger button will activate a pulse only when the device is in full contact with the skin.
- Verify that the treatment window is properly attached to your device (the device will not be operable if the treatment window is not attached or not properly attached to the device).
- Verify that the area you are trying to treat is not too small, curved or bony. The entire treatment window should touch the treatment area. Pulling the skin flat can help to achieve this.

- If the skin tone identified by the skin tone and contact sensor is too dark, the flash activation will be blocked, and the status indicator will turn "momentarily" red.

In case of general error:

- Reset the device by unplugging the power cord from the electrical outlet on the wall, wait 10 seconds and then reconnect it. Then try turning on the device after resetting.
- In case of dirt or debris on the treatment window, you should clean the treatment window according to the cleaning section. If you can't get the area clean, contact our customer service.
- In case of a pulse delay after repositioning the device on a different treatment area, you should consider performing the treatment in a cooler environment. The device may pulse at a slower rate in a warm environment.

If the problem is not solved, please contact our customer service.

Customer Service

For more information about the device, if your device is broken, stopped working, or if you need more assistance, please contact us:

 +44 0800-923-4833

 supportUK@sensica.com

Device specification parameters

Manufactured by:

ELG trade Ltd.

Emitted light spectrum:

Large treatment window:
(yellow optical filter)
475nm to 1,200nm max.

Precise treatment window:
(red optical filter)
550nm to 1,200nm max.

Pulse duration:

5 ± 2 [ms]

Emitted energy flux:

5 [J/cm²]

Selectable energy levels:

6

Treatment area:

Large treatment window:
4.5 cm²

Precise treatment window:
2 cm²

Speed:

1 pulse every 3 seconds

External Power Supply:

12[V], 1.8[A], 21.6[W]
100–240 VAC, 50/60Hz

Device Weight and Size:

16.3 [Oz], 7.9 X 6.9 X 2.7 [in]
461 [gr], 200 X 175 X 70 [mm]

Average active efficiency:

87.1%

Efficiency at low load (10%):

81.9%

Technology:

Home-use IPL

Package includes:

Hand-held device, Large treatment window, Precise treatment window, Power adapter and User guide.

Use the device only with the supplied external power and accessories. For your safety, never use DC adapters which are not approved by the manufacturer.

Storage Instructions:**Storage and Transport Conditions**

Temperature	-10°C to 70°C (14°F to 158°F)
Humidity	0 to 90% RH at 55°C (131°F)
Pressure	non-condensing

Operational Conditions

Temperature	10°C to 40°C (50°F to 104°F)
Humidity	30% to 70% RH
Max operating altitude	10,000 ft.

Labeling

Model: IPL300XXYYZZ

Double insulation

Consult Instructions for Use

Waste Electronic and Electrical Equipment

Keep away from eyes

Please read all warnings and safety precautions prior to use

Do not use the device on tanned skin or after sun exposure

Do not let the device become wet. Keep away from bathtubs, showers, basins, or other vessels containing water: Keep dry

Rechargeable battery

IP2X: Object of 12.5 mm diameter, shall not penetrate.

IP22

IPX2: Vertically falling drops shall have no harmful effects.

Warning: Optical radiation

Notes on handling batteries

If your skin or eyes come into contact with fluid from the battery cell, flush out the affected areas with water and seek medical assistance.

Protect batteries from excessive heat.

Risk of explosion! Never throw batteries into a fire.

Do not disassemble, open or crush the batteries.

Some warming of the device and power adapter during charging is normal.

Do not charge the battery in a room temperature above 40°C (104 F). The device is equipped with battery-over-heat protection and will not charge if the room temperature exceeds 40°C (104 F).

Never cover the appliance and the adapter during charging.

Only use chargers specified in the instructions for use.

Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in this user guide regarding correct charging must be observed at all times.

Fully charge the battery prior to initial use (see Device Set-up section above).

In order to achieve maximum battery service life, fully charge the battery at least twice a year.

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